

deSouffer
MEDICAL



VdriveTM

“BEST IN CLASS POWER TO WEIGHT RATIO”

THE ULTIMATE IN PERFORMANCE VETERINARY POWER TOOLS




THE Vdrive™ SYSTEM

Since the development of the first Slocum TPLO saw, in partnership with Dr Barclay Slocum DVM, the De Soutter Medical Veterinary range is now recognised as being the best in its class.

Our Vdrive™ veterinary range is lightweight and fully modular, and comprises a large and small bone system which share a common battery platform. The V-MBQ is suitable for large bone surgery, with a uniquely reliable TPLO attachment. Whilst the V-MBU is equally suitable for small bone surgery.

We are pleased to be introducing our elegant straight case TPLO saw with its unique dual position trigger function, enabling complete flexibility of grip.

A range of aseptic and sterile batteries with Li-Ion technology ensures a class-leading power to weight ratio.

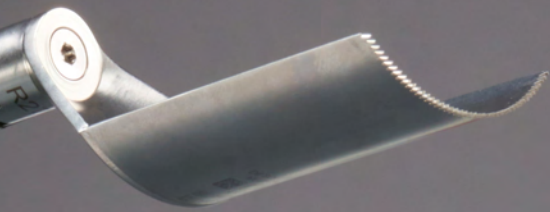

Vdrive™ V-MBU

SMALL BONE POWERED INSTRUMENT RANGE



V-MBU-470

Fitted with wire driver attachment

ADVANCED FUNCTION CONTROLS

- High speed and oscillating drilling.
- Low speed drilling, screw driving and tapping.
- Torque limiting for screw driving.
- Forward and reverse triggers.
- Mode selector button with mechanical safety lock.

POWER OPTIONS

- 14.6V Li-Ion aseptic batteries.
- 13.2V Li-Ion sterile batteries.
- Corded direct power.

V-MBU SYSTEM

HANDPIECE



ORDER CODE	MODEL	DESCRIPTION	SPEED (rpm)
1292444	V-MBU-470	Twin Trigger	0 - 1350

ATTACHMENTS

	ORDER CODE	MODEL	DESCRIPTION	SPEED (rpm)	CANNULATION Ø mm	
	18780	V-DU-470	0.5 - 4.0 mm (5/32") Keyed	0 - 1350	4.1	
	18770	V-DU-470	0.5 - 7.4 mm (1/4") Keyed	0 - 1350	4.2	
	18800	V-DU-470	A.O. Synthes/ASIF	0 - 1350	2.55	
	18790	V-DU-470	Hudson/Zimmer	0 - 1350	4.2	
	19480	V-BU-472	High speed bur Uses R Series burs	0 - 27000	-	
		Standard bur guard Order code 19490		Long bur guard Order code 19500		Extra long bur guard Order code 19510
	18840	V-RU-470	Hudson/Zimmer	0 - 350	4.2	
	18830	V-RU-470	A.O. Synthes	0 - 350	4.2	
	18820	V-RU-470	0.5 - 7.4 mm (1/4")	0 - 350	4.2	
	18870	V-SU-470	Mini Sagittal saw Uses S85 type blades	0 - 22000 cpm	-	
	18880	V-NU-470	Quick release sagittal saw Uses S88 type blades	0 - 20500 cpm	-	
	18890	V-CU-470	Small reciprocating saw Uses S92 type blades	0 - 8000 cpm	-	

	ORDER CODE	18910
	MODEL	V-WU-471
	DESCRIPTION	0.7 - 4.0 mm capacity
	SPEED (rpm)	0 - 1350
	CANNULATION Ø mm	4.0

	ORDER CODE	18900
	MODEL	V-WU-470
	DESCRIPTION	0.6 - 1.8 mm capacity
	SPEED (rpm)	0 - 1350
	CANNULATION Ø mm	1.8

	ORDER CODE	18920
	MODEL	V-PU-470
	DESCRIPTION	1.8 - 3.2 mm capacity
	SPEED (rpm)	0 - 1350
	CANNULATION Ø mm	3.2

Vdrive™ V-MBQ

LARGE BONE POWERED INSTRUMENT RANGE





V-MBQ-707

Fitted with slocum saw attachment and Jacobs drill.

ADVANCED FUNCTION CONTROLS

- A single handpiece will cover all large bone veterinary procedures with its comprehensive range of attachments.
- Mode selector button with mechanical safety lock.
- Both double and single trigger options available giving total flexibility.


HANDPIECE SPECIFICATION

	ORDER CODE	MODEL	DESCRIPTION	SPEED (rpm)
	1291774	V-MBQ-707	Single Trigger	0 - 1350
	1291884	V-MBQ-708	Twin Trigger	0 - 1350
* only to be used with De Soutter Medical V-MBQ attachments				


STRAIGHT CASE SYSTEM



DEDICATED HANDPIECE SPECIFICATION

	ORDER CODE	MODEL	DESCRIPTION	SPEED (cpm)
	1293884	V-DBO-800	Dedicated straight case TPLO saw Uses De Soutter Medical R series blade * Also accepts Synthes® blades Wrench order code 21150	0 - 13500

MODULAR HANDPIECE SPECIFICATION

	ORDER CODE	MODEL	DESCRIPTION	SPEED (rpm)
	1293774	V-MBQ-807*	Modular straight case handpiece * To be used with De Soutter Medical V-MBQ attachments	0 - 1350

V-MBQ ATTACHMENTS

	ORDER CODE	MODEL	DESCRIPTION	SPEED (rpm)	CANNULATION Ø MM
	18320	V-DQ-708	0.5 - 6.4 mm (1/4") keyed	0 - 1350	4.4
	18330	V-DQ-708	Small A.O Synthes/ASIF	0 - 1350	2.55
	18440	V-RQ-708	0.5 - 6.4 mm (1/4") keyed	0 - 330	4.4
	18290	V-KQ-707	Sagittal saw Uses S89 type blades	0 - 13500 cpm	
	18300	V-NQ-707	Small sagittal saw Uses S88 type blades	0 - 20000 cpm	
	18690	V-OQ-708*	TPLO oscillating saw Uses De Soutter Medical R series blade * Also accepts Synthes® blades Wrench order code 21150	0 - 13500 cpm	
	18270	V-OQ-708	TPLO oscillating saw Uses Slocum blades	0 - 13500 cpm	
	18310	V-WQ-707	0.7 - 4.0 mm wire capacity	0 - 1350	4.0

BATTERY AND ELECTRICAL COMPONENTS



Why choose Lithium Ion battery technology?

De Soutter Medical Lithium Ion batteries offer the surgical team and sterile services a number of tangible benefits over conventional battery technology, including:

- Li-Ion technology improves battery 'capacity to weight' ratio when compared to NiMH or NiCad chemistries.
- Significantly extends the usable run time of the power tool.
- Eliminates 'memory effect' normally associated with other battery types.
- Unique **fuel gauge** incorporated within De Soutter Medical Li-Ion batteries allows the surgical team to instantly assess the battery capacity - both before and during use.



STERILE BATTERIES - LI-ION

	ORDER CODE	MODEL	DESCRIPTION	VOLTAGE	CAPACITY (mAh)
	17210	SB-703	Small battery	13.2V	1100
	17220	SB-704	Large battery	13.2V	2500

ASEPTIC BATTERIES

SMALL BATTERY



Battery pack



Battery housing



Battery shield

ORDER CODE	MODEL	DESCRIPTION	VOLTAGE	CAPACITY (mAh)
17910	AB-703	Battery pack	10.95V	2000
18850	V-AH-703	Battery housing	-	-
18100	AS-703	Battery shield	-	-

LARGE BATTERY



Battery pack



Battery housing



Battery shield


ORDER CODE	MODEL	DESCRIPTION	VOLTAGE	CAPACITY (mAh)
17920	AB-704	Battery pack	14.6V	2000
18860	V-AH-704	Battery housing	-	-
18110	AS-704	Battery shield	-	-

BATTERY AND ELECTRICAL COMPONENTS

BC-700 MULTI BAY CHARGER AND MODULE OPTIONS

CHARGER ORDER CODE	CHARGER MODULE DESCRIPTION	CHARGER MODULE ORDER CODE	FOR BATTERY TYPE
15700	Without charger module fitted		
16230	With 4-off CA-703 charger modules	18120	AB-703 / AB-704
16220	With 4-off CA-700 charger modules	15750	SB-703 / SB-704

SINGLE BAY CHARGERS

	ORDER CODE	MODEL	FOR BATTERY TYPE
	17990	BC-706	SB-703/SB-704
	18000	BC-707	AB-703
	18010	BC-708	AB-704

MAINS CABLES

Mains power, 3 metres long

ORDER CODE	COUNTRY
8290	UK
8300	Euro
8310	US / Canada
8320	Australia / New Zealand
14700	Swiss

BC-700 BATTERY CHARGER

The BC-700 Battery Charger requires minimal user input, simply insert battery and the BC-700 will charge, recover and recondition it to its optimum state.



TPLO SAW BLADES

De Soutter Medical manufactured TPLO saw blades offer a range of advantages;

- Thinner blades - reduces weight, vibration and requires less cutting force
- Reduced cut thickness - minimises heat gain
- A sharper blade lasts longer
- Tooth offset improves cutting performance
- Tapered hub design reduces vibration and eliminates risk of blade coming undone
- Simple single captive screw makes fitting easier



BLADE RANGE / PART NO.

	Length 45mm					Length 50mm	
Radius	12	15	18	21	24	27	30
DSM/Synthes®	20160	20170	20180	20190	20200	20210	20220

Synthes® is a registered trademark of Synthes-Stratec. De Soutter Medical are not associated with the Synthes® trademark.

STERILISATION

	SYSTEM	Half size Container & insert tray Order code	Half size Insert tray Order code
	V-MBU ASEPTIC	19210	19190
	V-MBQ-70X ASEPTIC	19230	19220
	V-MBQ-807 ASEPTIC	20430	20380
	UNIVERSAL	19000	667923
	Sterile filtered container Accepts half size wire insert tray Order code: 19370		Filter paper - single use pack of 100 Order code: 657243
	Security tag, pack of 100 Order code: 657253		Sterile indicator label, pack of 100 Order code: 657263
	Lid gasket Order code: 669903		

CORDLESS BATTERY POWERED SAW

SAW SPECIFICATION



ORDER CODE	MODEL	VOLTAGE	OSCILLATING SPEED (CPM)	WEIGHT (GRAMS)
1286224	CSB-100	10.8 Vdc	16000	1050

Supplied with saw blade part no. 6250.

BATTERY AND BATTERY CHARGER

ORDER CODE	MODEL	VOLTAGE	SUPPLY FREQUENCY
15020	BC-100	100 - 240 Vac	50 - 60 Hz
15010		10.8 Vdc NiMH	

CSB-100-A CAST SAW SYSTEM



ORDER CODE	DESCRIPTION	ORDER CODE
1286774	CSB-100 Saw	1286224
	CB-100 Battery x 2	15010
	BC-100 Charger	15020
	Saw blade	6250
	Hexagon wrench	204903
	Carry case	15030

CSB-100

The De Soutter Medical CSB-100 saw allows the user total freedom of movement, which is particularly useful for the veterinary profession where the cordless concept is the only practical and affordable solution.



Vdrive™

“BEST IN CLASS POWER TO WEIGHT RATIO”

United Kingdom

(International Head Office)
De Soutter Medical Limited
Halton Brook Business Park,
Weston Road, Aston Clinton,
Aylesbury, Bucks, HP22 5WF
Tel: +44 (0) 1296 634 000
Email: info@de-soutter.com

Belgium

De Soutter Medical Belgium
Bessemmerstraat 14, 3620 Lanaken
Tel: +32 (0) 89/47 15 37
Email: belgium@de-soutter.com

Italy

De Soutter Medical Italia
Località Fornace SNC
27022 Casorate Primo - PV
Tel: +39 (0) 2 9009 4098
Email: italy@de-soutter.com

Australia

De Soutter Medical Australia Pty Ltd
2/12-14 Apollo Drive, Hallam,
Victoria 3803
Tel: +61 (0) 3 9702 4441
Email: australia@de-soutter.com

France

De Soutter Médical France
1 252 Avenue du Parc des Expositions
33260 La Teste de Buch
Tél. : +33 (0) 5 56 54 89 36
Email : france@de-soutter.com

Netherlands

De Soutter Medical Netherlands
Gelderlandhaven 2X
3433 PG Nieuwegein
Tel: +31 (0) 85 0491480
Email: nederland@de-soutter.com

Austria

De Soutter Medical Austria
Zweigniederlassung Österreich
Hauptstrasse 122-124/Büro 01
3420 Kritzensdorf
Tel: +43 (0) 2243 239000
Email: austria@de-soutter.com

Germany

De Soutter Medical Germany
Bahnhofstraße 4, Kladenfloss
D-66625 Nohfelden
Tel: +49 (0) 68 52-99 12 46
Email: deutschland@de-soutter.com

United States of America

De Soutter Medical USA, Inc.
224 Rolling Hill Road, Suite 12A
 Mooresville, NC 28117
Tel: +1 (704) 655 9040
Email: usa@de-soutter.com

211-EN-12

deSoutter
MEDICAL