

deSoutter
MEDICAL



sternudriveTM

www.de-soutter.com

The De Soutter Medical MBC sternudrive Lite™ is a powerful, yet versatile modular sternum saw system for primary, revision and paediatric surgery. The handpiece is powered by either sterile or aseptic Li-ion batteries, assuring class leading power to weight ratio.



MBC-709
sternudrive Lite™
Handpiece



Oscillating
sagittal
sternum saw
Uses S87 saw blades



Oscillating
paediatric sagittal
sternum saw
Uses S88-9xx saw blades



Reciprocating
sternum saw
Uses S84-4xx saw blades



Sternum guard



Oscillating sternum saw
Uses S88-9xx
saw blades

STRAIGHT CASE MODULAR STERNUM SAW

The De Soutter Medical MBC-809 & DBR-800 sternudrive Lite™ are elegant straight case handpieces with a unique dual position trigger function, allowing complete grip flexibility.



MBC-809
sternudrive Lite™
Straightcase Handpiece

	Model	Order code	Description	Speed cpm/rpm
	MBC-709	1291554	Modular handpiece Cannulation 4.4 mm	-
	MBC-809	1293664	Modular straight case handpiece	-
	CQ-709	18380 16490	Reciprocating sternum saw Sternum guard	0 - 13000
	NQ-709	18450	Oscillating paediatric sternum saw Uses S88-9xx saw blades	0 - 13000
	KQ-709	17650	Oscillating sagittal sternum saw Uses S87 saw blades	0 - 13000
	DQ-709	19380	0.5 - 4.0 mm (5/32") keyed drill Cannulation 4.4 mm	0 - 1300
	DQ-709	19430	Small A.O. Synthes/ASIF drill Cannulation 2.1 mm	0 - 1300
	OQ-709	21390	Oscillating sternum saw Uses S88-9xx saw blades	0 - 13000
	WQ-709	19390	Wire driver 0.7 - 4.0 mm	-
	Model	Order code	Description	Speed cpm/rpm
	DBR-800	1295334	Oscillating Saw straight case handpiece Uses S88 saw blades	0 - 13500

DB BATTERY STERNUM SAWS

The DBC-470 has a unique active damping system which reduces vibration by 30% and noise by up to 40%. A single piece micro adjustable blade guard provides optimum control of the cutting direction for primary sternotomies.

DBR-700 for revision sternotomies

DBR-470 paediatric

DBC-470 for primary sternotomies

Simple twist and click sternum blade guard for quick and efficient use

Sternum guard
Order code 16490

Model	Order code	Speed cpm	Stroke
DBC-470	1290224	0 - 14000	3.7 mm
DBR-470	1291664	0 - 15000	-
DBR-700	1290884	0 - 13500	-

DP AIR STERNUM SAWS

DPC-470 offers the same innovative functionality as the battery sternum saw. A single piece, twist and click, micro adjustable blade guard provides quick and efficient use with optimum control of the cutting direction. For revision sternotomies, the DPX-170 oscillating saw handpiece is also offered.



DPC-470
for primary
sternotomies

Wash assist air hose



This hose system is recommended where machine washing of air hoses is required. By releasing the probe end from the hose, which is retained on the hose via the rubber strap, the two hose ends can be joined and sealed prior to machine washing.

Probe end	Model	3 m hose	4 m hose	5 m hose
Order code				
MA7 Schrader	DPC-470	13300	13400	13500
	DPX-170	13000	13100	13200
AO	DPC-470	13310	13410	13510
	DPX-170	13010	13110	13210
Draeger	DPC-470	13330	13430	13530
	DPX-170	13030	13130	13230
AGA Inst.	DPC-470	13320	13420	13520
	DPX-170	13020	13120	13220

Other probe connectors are available on request.

Air inlet washing cap.



Order code 9950 for DPC-470
Order code 608073 for DPX-170



Model	Order code	Speed cpm	Stroke	Operating pressure	Sternum guard
DPC-470	226614	0 - 16000	3,7 mm	6 - 7 bar	Order code 16490
Oscillating arc					
DPX-170	223224	0 - 17000	6°	6 - 7 bar	



BATTERY CHARGERS



BC-700 BATTERY CHARGER



shown with aseptic batteries.

The BC-700 Battery Charger requires minimal user input, simply insert battery and the BC-700 will charge, recover and recondition it to its optimum state. **The large touch-screen display** with intuitive menu navigation, displays instant battery charge status and condition. Includes a multi-lingual display and other user settings.

The BC-700 also offers a **unique healthcheck facility**, profiling the condition of an individual battery, automatically reconditioning it, to provide optimal capacity. Alternatively, a manual healthcheck can be easily selected. The BC-700 Battery Charger accepts all De Soutter Medical batteries, ensuring total backwards compatibility.



BC- 700 MULTI BAY CHARGER AND MODULE OPTIONS

Charger - Order code	Charger module description	Charger module Order code	For battery type
15700	Without charger module fitted	-	-
16230	With 4-off CA-703 charger modules	18120	 AB-703/AB-704
16220	With 4-off CA-700 charger modules	15750	 SB-703/SB-704

Mains power cordsets - 3 metres long.
Please specify - Order code

UK	8290
Euro	8300
US/Canada	8310
Australia/New Zealand	8320
Swiss	14700
Brazil	21340

The BC-700 charger modules can be quickly and easily changed to suit the batteries being charged.



SINGLE BAY CHARGERS

Model	Battery type - Order code	For battery type
BC-706	17990	SB-703/SB-704
BC-707	18000	AB-703
BC-708	18010	AB-704

LITHIUM ION STERILE BATTERIES

Why choose Lithium Ion battery technology



De Soutter Medical Lithium Ion batteries offer the surgical team and sterile services a number of tangible benefits over conventional battery technology.

Key advantages are:

- Li-Ion technology improves battery 'capacity to weight' ratio when compared to NiMH or NiCad chemistries.
- Significantly extends the useable run time of the power tool.
- Eliminates 'memory effect' normally associated with other battery types.
- De Soutter Medical sterile Li-Ion technology readily withstands repeated machine washing and autoclaving.
- Unique fuel gauge incorporated within De Soutter Medical Li-Ion batteries allows the surgical team to instantly assess the battery capacity - both before and during use.



LITHIUM ION STERILE BATTERIES

	Model	Description	Order code	Voltage	Capacity mAh
	SB-703	Small battery	17210	13.2V	1100
	SB-704	Large battery	17220	13.2V	2500



SMALL BATTERY

	Battery pack		Aseptic housing		Aseptic shield
Model	Description	Order code	Voltage	Capacity mAh	
AB-703	Battery pack	17910	10,95V	2000	
AH-703	Battery housing	18080	-	-	
AS-703	Battery shield	18100	-	-	

LARGE BATTERY

	Battery pack		Aseptic housing		Aseptic shield
Model	Description	Order code	Voltage	Capacity mAh	
AB-704	Battery pack	17920	14,6V	2000	
AH-704	Battery housing	18090	-	-	
AS-704	Battery shield	18110	-	-	

LOADING AN ASEPTIC BATTERY PACK



1

Press both catches to open housing.




2

With the aseptically housing open, the sterile operator fits the aseptically shield over the opening.



3

The non-sterile operator inserts the fully charged aseptically battery into the housing.

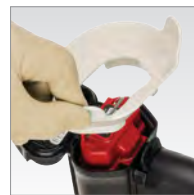

4

The non-sterile operator removes the aseptically shield.



5

The sterile operator closes the aseptically housing until the lid clicks shut, ensuring the housing is fully closed.


SAW BLADES

RECIPROCATING BLADES FOR DBC-470, DPC-470 AND CQ-709




Blade Order code	Blade thickness mm	Cut thickness mm	Tooth pitch mm	
S84-400	0.80	1.2	1.4 (18 TPI)	Uni directional cut
S84-401	0.80	1.2	1.4 (18 TPI)	Bi directional cut
S84-402	0.60	1.2	1.4 (18 TPI)	Uni directional cut
S84-406	0.80	1.0	1.3	

OSCILLATING BLADES FOR OI-280 - FOR USE WITH MCI-28X




Blade Order code	Blade thickness mm	Cut thickness mm	Tooth pitch mm	Useable length mm	Blade width mm
S86-5695	0.38	0.6	1.4	31	50

OSCILLATING BLADES FOR DBR-700 AND KQ-709



Blade Order code	Blade thickness mm	Cut thickness mm	Tooth pitch mm	Useable length mm	Blade width mm
S87-112	0.6	0.8	1.9	65	50
S87-113	0.6	0.8	1.9	65	50
S87-118 Paediatric	0.51	0.9	1.4	23	50
S87-119 Paediatric	0.51	0.9	1.4	23	63
S87-120	0.51	1.0	1.6	35	82
S87-121	0.51	1.0	1.6	35	82
S87-123	0.51	0.9	1.4	23	50
S87-125	0.6	0.8	1.2	48	43
S87-512	0.6	0.6	2.1	65	50

OSCILLATING BLADES FOR DBR-470, DBR-800, NQ-709 AND OQ-709



Blade Order code	Blade thickness mm	Cut thickness mm	Tooth pitch mm	Useable length mm	Blade width mm
S88-901	0.51	0.8	1.6	37	25
S88-905	0.51	0.8	1.6	19	13 & 62
S88-908	0.60	0.6	1.4	19	50
S88-909	0.60	0.6	1.4	19	63
S88-918	0.60	0.6	1.4	25	50
S88-919	0.60	0.6	1.4	25	63
S88-920	0.51	0.8	1.6	37	82

For use with DBR-800.

STERILISATION INSERT TRAYS

The introduction of a high performance hard anodising technology onto our handpieces has provided:

- Excellent corrosion resistance and protection against highly alkaline detergents. Durability and long term reliability is significantly improved.
- Ultra tough and exceptionally scratch resistant.



ASEPTIC BATTERY STERILISATION SYSTEM

Set up - optional tray contents	Container and insert tray Order code	Insert tray Order code
1 off DBR-470 or DBC-470 1 off sternum guard (DBR_470) 1 off aseptic housing and shield	18210	18200
1 off DBR-700 handpiece 1 off aseptic housing and shield	19170	19070
1 off DBR-800 handpiece 1 off aseptic housing and shield	21370	21350
1 off MBC-709 handpiece 3 off attachments 1 off sternum guard 1 off aseptic housing and shield	19280	19300
1 off MBC-809 handpiece 1 off MBC-709 handpiece 6 off attachments 1 off sternum guard 1 off aseptic housing and shield	20410	20360

STERILE BATTERY STERILISATION SYSTEM

Set up - optional tray contents	Container and insert tray Order code	Insert tray Order code
1 off DBR-470 or DBC-470 1 off sternum guard (DBR_470) 1 off sterile battery	19270	19110
1 off DBR-700 handpiece 1 off sterile battery	17890	17880
1 off DBR-800 handpiece 1 off sterile battery	21380	21360
1 off MBC-709 handpiece 3 off attachments 1 off sternum guard 1 off sterile battery	19200	19100
1 off MBC-809 handpiece 1 off MBC-709 handpiece 6 off attachments 1 off sternum guard 2 off sterile batteries	20420	20370
1 off DPC-470 handpiece 1 off sternum guard	17160	17170
2 off SB-70X batteries	17390	17400

Wire container for aseptic and sterile trays
Order code 19000

STERILISATION CONTAINERS

WIRE CONTAINERS



Standard wire container
Order code 19010
Accepts standard insert tray

Half size wire container
Order code 19000
Accepts half size insert tray

STERILE FILTERED CONTAINERS

Accepts standard insert tray



Standard sterile filtered container
Order code 19250

PTFE sterile filtered container
Order code 19810

Lid gasket
Order code 670023

PTFE filter cassette
Order code 671073



Accepts half size insert tray
Half size sterile filtered container
Order code 19370
Lid gasket
Order code 669903



Filter paper - single use
pack of 100
Order code 657243



Security tag
pack of 100
Order code 657253



Sterile indicator label, pack of 100. Order code 657263

UNITED KINGDOM (INTERNATIONAL HEAD OFFICE)

De Soutter Medical Limited
Halton Brook Business Park, Weston Road,
Aston Clinton, Aylesbury, Bucks HP22 5WF
Tel: +44 (0) 1296 634 000
Email: info@de-soutter.com

AUSTRALIA

De Soutter Medical Australia Pty Ltd
2/12-14 Apollo Drive, Hallam
Victoria 3803
Tel: +61 (0) 3 9702 4441
Fax: +61 (0) 3 9702 4484
Email: australia@de-soutter.com

AUSTRIA

De Soutter Medical Austria
Hauptstrasse 122-124/Büro 01
3420 Kritzensdorf
Österreich
Tel: +43 (0) 2243 239000
Email: austria@de-soutter.com

BELGIUM

De Soutter Medical Belgium
Bessemersstraat 14
3620 Lanaken
Tel: +32 (0) 89/47 15 37
Email: belgium@de-soutter.com

FRANCE

De Soutter Médical France
1 252 Avenue du Parc des Expositions
33260 La Teste de Buch
Tél. : +33 (0) 5 56 54 89 36
Email : france@de-soutter.com

GERMANY

De Soutter Medical Deutschland
Bahnhofstraße 4
D-66625 Nohfelden
Tel: +49 (0) 68 52-99 12 46
Email: deutschland@de-soutter.com

ITALY

De Soutter Medical Italia
Località Fornace SNC
27022 Casorate Primo - PV
Tel: +39 (0) 2 9009 4098
Email: italy@de-soutter.com

NETHERLANDS

De Soutter Medical Nederland
Gelderlandhaven 2X
3433 PG Nieuwegein
Tel: +31 (0) 85 0491480
Email: nederland@de-soutter.com

UNITED STATES OF AMERICA

De Soutter Medical USA, Inc.
224 Rolling Hill Road, Suite 12A
 Mooresville, NC 28117
Tel: +1 (704) 655 9040
Email: usa@de-soutter.com

deSoutter
MEDICAL

For further information on De Soutter Medical powered instruments and cast saw systems please contact the address in your area.